

clinell®

Drain Disinfectant

Proven protection against high
risk and hard-to-kill organisms





Product introduction	4
Wet surface biofilms	6
Current techniques aren't working	9
Peracetic acid advantage	10
Order info	11

Contaminated wet and dry surfaces contribute to the transmission of pathogens associated with healthcare-associated infections (HCAI)^{1,2}.

Over half of Intensive Care Units (ICU) have a sink(s) contaminated by multidrug resistant bacteria.*

In clinical settings, sinks meant for hand hygiene are used for waste disposal (including clinical fluid and drinks); hand washing accounted for only 4% of uses in one study⁸.

Contaminated sinks and drains are an important factor in the transmission of key Gram-negative bacteria, including *Pseudomonas aeruginosa* and carbapenemase-producing *Enterobacteriaceae* (CPE)^{3,4}.

Two thirds of ICUs have a sink(s) that are contaminated with a CPE or an IRPA.*†

6 in 10 ICU sinks are inadequately disinfected or not disinfected at all.*



*Based on a multicentre study of 73 French ICUs wards covering 996 beds²⁹.
† Imipenem-resistant *Pseudomonas aeruginosa*

High-performance disinfection

Contamination plays an important role in the transmission of HCAs¹. It is common for sinks and drains to be contaminated with antibiotic resistant bacteria⁶.

Some of the most persistent outbreaks come from drains. In these ideal conditions, microorganisms form biofilms – protective structures rendering traditional disinfectants ineffective.

Clinell Drain Disinfectant is a high-performance disinfectant that provides proven protection against high-risk and hard-to-kill organisms.

More effective than traditional disinfectants

Getting to the source of the problem, Clinell Drain Disinfectant wipes out biofilms throughout the drainage system.

Simple to use oxidative technology


One sachet of Clinell Drain Disinfectant granules, combined with water, produces enough powerful oxidative agents including peracetic acid (see page 10) to break down biofilms and kill the microbes living within them.

Strong on microbes, gentle on drains

Despite generating strong oxidising agents, our formula maintains a near neutral pH. Regular use will prevent biofilm regrowth while conserving the integrity of the drain, preventing outbreaks before they happen.

High-risk & hard-to-kill organisms

Different microorganisms have different physiological structures and therefore exhibit differing tolerance to external factors such as chemical disinfection. Clinell Drain Disinfectant is proven to provide a high level of efficacy against otherwise hard-to-kill microbes including bacteria, viruses, fungi and spores.

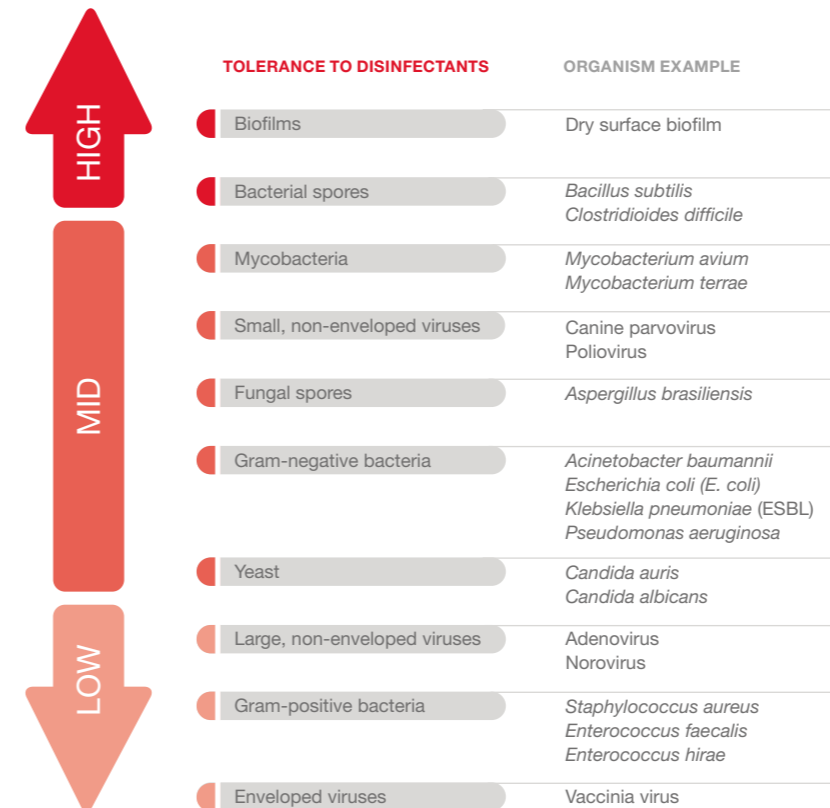

Typical tolerance of microorganism types to disinfectants, adapted from McDonnell & Russell²⁵

Wet surface biofilms

Sink and shower drains provide an ideal environment for microbes to form biofilms – cities of microbes living within a protective layer. Biofilms provide increased resistance to traditional disinfectants²⁶⁻²⁸, indicating why some outbreaks may seem impossible to resolve.

Even in the absence of outbreaks and known cases of infection, sinks have been found to be contaminated with microbes that have the potential for causing an outbreak⁷. Activities including hand hygiene supply bacteria that colonise the drainage system, and disposed fluids provide nutrients, thus supporting growth of biofilms⁹.


Biofilms form when free-floating microbes attach to a surface and change state. The secretion of Extracellular Polymeric Substances (EPS) creates a protective 'biofilm matrix' that traditional disinfectants such as chlorine can't penetrate.

Protected by the matrix, the microbes living within the protective layer are sheltered from external threats and are able to 'exchange genes' and transfer genetic material. This transfer of genetic material can occur between the dozens of pathogenic species that live within the biofilm. Resistance in one species can then transfer and can give rise to new antibiotic resistant bacteria^{10,11}.

Planktonic ("free-swimming") microbes attach to a surface to begin forming a biofilm.

The microbes produce a mix of Extracellular Polymeric Substances (EPS) – the "biofilm matrix" providing an added layer of microbial defence from disinfectants.


Inside the biofilm microorganisms are free to transfer antibiotic resistance genes.

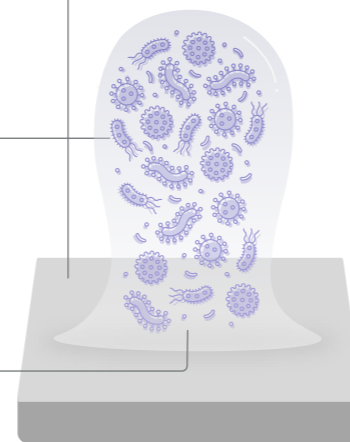

Contaminating the clinical setting

Once established, the wet surface biofilm continues to grow. When the area is treated with traditional disinfectants, the free-floating microbes outside of the biofilm are quickly killed. Microbiological samples taken will show the surface to be clean and disinfected. However, this is not the case. Inside the biofilm, the microbes survive and increase in numbers. Eventually, the biofilm will seed free-floating microbes back into the environment.

Laboratory experiments have shown that contamination in sinks and drains can potentially be transferred to the hands of healthcare workers and, subsequently, to patients^{12,13}. These papers have been followed by clinical studies demonstrating exactly that occurrence¹⁴⁻¹⁷.

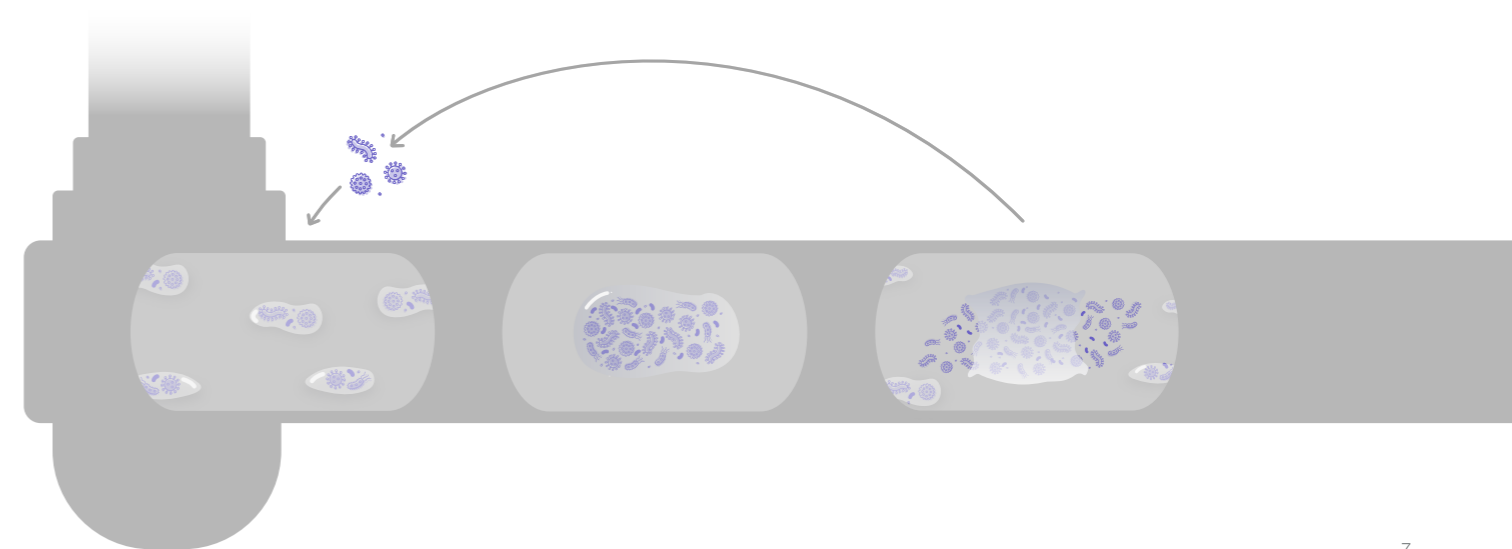

Spread from the sink to the patient: In situ study using green fluorescent protein (GFP)-expressing *Escherichia coli* to model bacterial dispersion from hand-washing sink-trap reservoirs¹³.

Kotay et al. *Applied and Environmental Microbiology*. 2017;83:e03327-16.

A hand washing sink model was established in a lab setting. The dispersal of green fluorescent protein (GFP)-expressing *Escherichia coli* bacteria was measured under various conditions. A regular supply of nutrients (to model typical use of sinks in healthcare settings) resulted in the formation of biofilm. Once the *E. coli* biofilm had developed, the GFP-expressing *E. coli* was dispersed widely around the sink. This model illustrates how biofilm-related contamination in a sink drain can be spread in a clinical setting.

Wet surface biofilms and healthcare-associated infections

Research has established that biofilms existing within drainage systems are linked to outbreaks and the spread of contamination in healthcare settings.

Wastewater drains: epidemiology and interventions in 23 carbapenem-resistant organism outbreaks⁵.

Carling PC. *Infection Control & Hospital Epidemiology*. 2018 39: 972-979.

This review of 23 outbreaks of carbapenem-resistant organisms summarises the recent evidence that contaminated drains and wastewater play a role in the continuation of these outbreaks.

Wastewater and drain-associated outbreaks were characterised by:

- Low density of new cases with long time periods between them
- Challenges with outbreak detection and definition
- Frequent colonisation of sinks and drains
- Apparent transfer of genetic material in drain biofilm
- Need for frequent drain disinfection in order to effectively tackle drain-related transmission

In summary, new genetic tools combined with new insights into the microbial ecology of biofilms provides evidence that contaminated drains and wastewater are an important factor in the transmission of carbapenem-resistant organisms in healthcare settings.



A prospective multicenter surveillance study to investigate the risk associated with contaminated sinks in the intensive care unit²⁹.

Valentin et al. *Clinical Microbiology and Infection*. (In Press) DOI: 0.1016/j.cmi.2021.02.018

The study aimed to assess the incidence of sink contamination by multidrug-resistant (MDR) *Pseudomonas aeruginosa* and *Enterobacteriaceae*, risk factors for sink contamination and splashing, and their association with clinical infections in the intensive care setting. From the 73 ICUs participating in the study, 50.9% (606/1191) of all sinks were contaminated by MDR bacteria.

These included:

- 41.0% of sinks used only for handwashing
- 55.3% of those used for waste disposal
- 23.0% of sinks treated daily with chlorine
- 62.0% of those untreated

A total of 459 sinks showed visible splashes and 30.5% were close to the bed (<2 m) with no barrier around the sink making them susceptible to splashing. MDR-associated bloodstream infection incidence rates were also examined.

The authors concluded that there were frequent and multifactorial infectious risks associated with contaminated sinks in ICUs.

Current techniques aren't working

Biofilms are hardy. Microbes living inside the biofilm can survive traditional disinfection, recover and seed microbes back into the environment.

Several different approaches have been taken to tackle contamination of sinks and drains in clinical settings. Changes in the structure of the clinical environment to reduce contact between staff/patients and sinks/drains, physical methods to remove biofilm and enhanced chemical disinfection¹⁸⁻²¹.

Techniques employing traditional disinfectant eradicate free-floating microbes in drainage systems but despite diligent disinfection, outbreaks still occur.

Studies in the clinical environment show that improved management of sinks and drains results in significantly less transmission of CPE¹⁸.

Intensive care unit wastewater interventions to prevent transmission of multispecies *Klebsiella pneumoniae* carbapenemase-producing organisms¹⁸.

Mathers et al. *Clinical Infectious Diseases*. 2018;67:171-178.

The study evaluated the impact of introducing covers for hoppers (a toilet-like waste disposal unit) and sink trap heater/vibration devices to heat-sanitise the drain and dislodge the formation of biofilm in an ICU. Decreases for all CPE acquisitions (a 49% reduction; P = .003) and CPE-positive clinical cultures (a 71% reduction; P < .001) per admission in patients exposed to an intervention unit were observed.

This suggests that improvement management of sinks and drains improved patient clinical outcomes.



Peracetic Acid advantage

Peroxygen chemicals such as peracetic acid are good candidates for a drain disinfectant because they are not prone to break down by dirt or organic matter, are safe for staff to use and have rapid biocidal activity (including sporicidal activity)²².

Peracetic acid has the noteworthy property of targeting both microbes in biofilms and the biofilm matrix itself – important in biofilm-rich environments such as sinks and drains⁹.

Clinell Drain Disinfectant combines with water to generate peracetic acid. It provides an effective option to tackle contaminated sink and shower drains, reduce risk of transmission and maximise patient safety.



It's a trap! The development of a versatile drain biofilm model and its susceptibility to disinfection²⁴.

Ledwoch et al. *Journal of Hospital Infection*. 106(4): 757-764.

Researchers at Cardiff University developed a novel in-vitro biofilm model to address the need for a robust, reproduceable and simple testing methodology for disinfection efficacy against a complex drain biofilm.

The efficacy of sodium hypochlorite 1000 ppm (NaOCl), sodium dichloroisocyanurate 1000 ppm (NaDCC), non-ionic surfactant (NIS) and peracetic acid 4000 ppm (PAA) was explored, simulating normal sink usage conditions.

Bacterial viability and recovery following a series of 15-min treatments were measured in three distinct parts of the drain. The drain biofilm consisted of 119 mixed species of Gram-positive and -negative bacteria.

NaOCl produced a >4 log₁₀ reduction in viability in the drain front section alone, NIS and NaDCC failed to control the biofilm in any drain sections.

Only the PAA formulation was able to significantly affect (>4 log₁₀ reduction) the biofilm throughout the drain model and importantly, prevent biofilm regrowth for a minimum of four days.

A systematic evaluation of a peracetic-acid-based high performance disinfectant²².

Humphreys et al. *Journal of Infect Prevention*. 2013;14:126-131.

This study evaluated the potential for peracetic-acid-based (PAA) disinfectants (such as Clinell Drain Disinfectant) to act as high performance disinfectants in healthcare settings.

When tested against bacteria and spores, PAA provided similar or improved performance than chlorine, especially when organic challenge was present or when tackling dried surface contamination.

These results suggest that PAA generating products provide an improved alternative to chlorine-based products.

Order info



Drain Disinfectant

24 sachets per box

Product code: CSDD24
NHSSC code: MKD85041



Indicator Tape

Single unit

Product code: CSDDT90
NHSSC code: MKD85040

To find out more, speak to your GAMA Healthcare Area Manager or visit www.gamahealthcare.com

References

1. Otter JA, Yezli S, French GL. The role played by contaminated surfaces in the transmission of nosocomial pathogens. *Infect Control Hosp Epidemiol*. Jul 2011;32(7):687-699.
2. Mitchell BG, Dancer SJ, Anderson M, Dehn E. Risk of organism acquisition from prior room occupants: a systematic review and meta-analysis. *J Hosp Infect*. Nov 2015;91(3):211-217.
3. Weingarten RA, Johnson RC, Conlan S, et al. Genomic Analysis of Hospital Plumbing Reveals Diverse Reservoir of Bacterial Plasmids Conferring Carbapenem Resistance. *mBio*. Feb 6 2018;9(1).
4. Berrouane YF, McNutt LA, Buschelman BJ, et al. Outbreak of severe *Pseudomonas aeruginosa* infections caused by a contaminated drain in a whirlpool bathtub. *Clin Infect Dis*. Dec 2000;31(6):1331-1337.
5. Carling PC. Wastewater drains: epidemiology and interventions in 23 carbapenem-resistant organism outbreaks. *Infect Control Hosp Epidemiol*. Aug 2018;39(8):972-979.
6. Buchan BW, Graham MB, Lindmair-Snell J, et al. The relevance of sink proximity to toilets on the detection of *Klebsiella pneumoniae* carbapenemase inside sink drains. *Am J Infect Control*. Jan 2019;47(1):98-100.
7. Lemarie C, Legeay C, Kouatchet A, et al. High prevalence of contamination of sink drains with carbapenemase-producing Enterobacteriaceae in 4 intensive care units apart from any epidemic context. *Am J Infect Control*. Feb 2020;48(2):230-232.
8. Grabowski M, Lobo JM, Gunnell B, et al. Characterizations of handwashing sink activities in a single hospital medical intensive care unit. *J Hosp Infect*. Nov 2018;100(3):e115-e122.
9. Kotay SM, Parikh HJ, Barry K, et al. Nutrients influence the dynamics of *Klebsiella pneumoniae* carbapenemase producing enterobacteriales in transplanted hospital sinks. *Water Res*. Jun 1 2020;176:115707.
10. Constantinides B, Chau KK, Quan TP, et al. Genomic surveillance of *Escherichia coli* and *Klebsiella* spp. in hospital sink drains and patients. *Microb Genom*. Jul 2020;6(7).
11. Muzslay M, Moore G, Alhussaini N, Wilson AP. ESBL-producing Gram-negative organisms in the healthcare environment as a source of genetic material for resistance in human infections. *J Hosp Infect*. Jan 2017;95(1):59-64.
12. Aranega-Bou P, George RP, Verlander NQ, et al. Carbapenem-resistant Enterobacteriaceae dispersal from sinks is linked to drain position and drainage rates in a laboratory model system. *J Hosp Infect*. May 2019;102(1):63-69.
13. Kotay S, Chai W, Guilford W, Barry K, Mathers AJ. Spread from the Sink to the Patient: In Situ Study Using Green Fluorescent Protein (GFP)-Expressing *Escherichia coli* To Model Bacterial Dispersion from Hand-Washing Sink-Trap Reservoirs. *Appl Environ Microbiol*. Apr 15 2017;83(8).
14. Jung J, Choi HS, Lee JY, et al. Outbreak of carbapenemase-producing Enterobacteriaceae associated with a contaminated water dispenser and sink drains in the cardiology units of a Korean hospital. *J Hosp Infect*. Apr 2020;104(4):476-483.
15. Heireman L, Hamerlinck H, Vandendriessche S, et al. Toilet drain water as a potential source of hospital room-to-room transmission of carbapenemase-producing *Klebsiella pneumoniae*. *J Hosp Infect*. Oct 2020;106(2):232-239.
16. Kotay SM, Donlan RM, Ganim C, Barry K, Christensen BE, Mathers AJ. Droplet- Rather than Aerosol-Mediated Dispersion Is the Primary Mechanism of Bacterial Transmission from Contaminated Hand-Washing Sink Traps. *Appl Environ Microbiol*. Jan 15 2019;85(2).
17. Hajar Z, Mana TSC, Cadnum JL, Donskey CJ. Dispersal of gram-negative bacilli from contaminated sink drains to cover gowns and hands during hand washing. *Infect Control Hosp Epidemiol*. Apr 2019;40(4):460-462.
18. Mathers AJ, Vegesana K, German Mesner I, et al. Intensive Care Unit Wastewater Interventions to Prevent Transmission of Multispecies *Klebsiella pneumoniae* Carbapenemase-Producing Organisms. *Clin Infect Dis*. Jul 2 2018;67(2):171-178.
19. Smolders D, Hendriks B, Rogiers P, Mul M, Gordts B. Acetic acid as a decontamination method for ICU sink drains colonized by carbapenemase-producing Enterobacteriaceae and its effect on CPE infections. *J Hosp Infect*. May 2019;102(1):82-88.
20. de Jonge E, de Boer MGJ, van Essen EHR, Dogterom-Ballering HCM, Veldkamp KE. Effects of a disinfection device on colonization of sink drains and patients during a prolonged outbreak of multidrug-resistant *Pseudomonas aeruginosa* in an intensive care unit. *J Hosp Infect*. May 2019;102(1):70-74.
21. Aranega-Bou P, Cornbill C, Verlander NQ, Moore G. A splash-reducing clinical handwash basin reduces droplet-mediated dispersal from a sink contaminated with Gram-negative bacteria in a laboratory model system. *J Hosp Infect*. Apr 22 2021.
22. Humphreys PN, Finan P, Rout S, et al. A systematic evaluation of a peracetic-acid-based high performance disinfectant. *J Infect Prev*. 2013;14(4):126-131.
23. Jones LD, Mana TSC, Cadnum JL, et al. Effectiveness of foam disinfectants in reducing sink-drain gram-negative bacterial colonization. *Infect Control Hosp Epidemiol*. Dec 5 2019:1-6.
24. Ledwoch K, Robertson A, Lauran J, Norville P, Maillard JY. It's a trap! The development of a versatile drain biofilm model and its susceptibility to disinfection. *J Hosp Infect*. Dec 2020;106(4):757-764.
25. McDonnell G, Russell AD. Antiseptics and disinfectants: Activity, action, and resistance. *Clin Microbiol Rev*. 1999;12(1):147-179. doi:10.1128/cmr.12.1.147
26. Condell O, Iversen C, Cooney S, et al. Efficacy of biocides used in the modern food industry to control *Salmonella enterica*, and links between biocide tolerance and resistance to clinically relevant antimicrobial compounds. *Appl Environ Microbiol*. 2012;78(9):3087-3097. doi:10.1128/AEM.07534-11
27. Leung CY, Chan YC, Samaranyake LP, Seneviratne CJ. Biocide resistance of *Candida* and *Escherichia coli* biofilms is associated with higher antioxidant capacities. *J Hosp Infect*. 2012;81(2):79-86. doi:10.1016/j.jhin.2011.09.014
28. Smith K, Hunter IS. Efficacy of common hospital biocides with biofilms of multi-drug resistant clinical isolates. *J Med Microbiol*. 2008;57(8):966-973. doi:10.1099/jmm.0.47668-0
29. Valentin, A. S., S. D. Santos, F. Goube, R. Gimenes, M. Decalonne, L. Mereghetti, C. Daniau, N. van der Mee-Marquet and S. I. group (2021). "A prospective multicentre surveillance study to investigate the risk associated with contaminated sinks in the intensive care unit." *Clin Microbiol Infect*. (In Press) DOI: 0.1016/j.cmi.2021.02.018

Use biocides safely. Always read the label and product information before use.